


InterGuide Tracheal Tube Introducers (Bougies) and InterForm Intubating Stylets

For use in the Difficult Airway

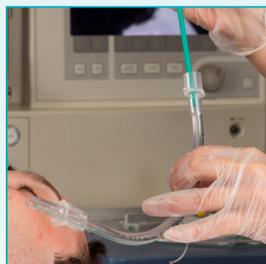

Introducing InterGuide and InterForm


The InterGuide and InterForm ranges have been designed to complement the InTube range of endotracheal tubes. Both product ranges are single use, latex free and supplied sterile in easy to open individual packs.

Features and Benefits

InterGuide

InterGuide is a flexible tracheal tube introducer commonly known as a bougie. It allows positive location of the trachea and subsequent placement of the endotracheal tube in difficult airway situations.

Atraumatic coude tip to reduce the potential for patient trauma

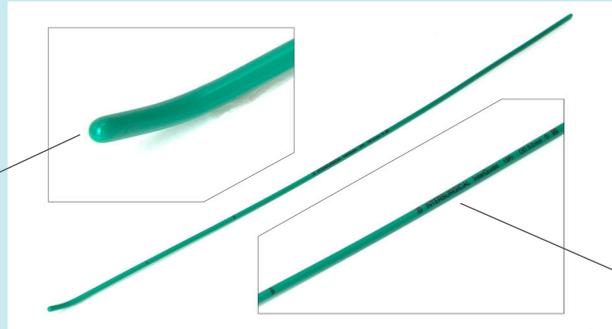

Low-friction surface for easy insertion and operation

Available in three sizes: 6FR, 10FR and 15FR

InterForm

InterForm is a malleable endotracheal stylet, allowing the user to form the endotracheal tube into a suitable shape to ease insertion and give more control of the airway device.

Aluminium core (MRI safe) allows user to easily form suitable shapes



Low-friction surface for easy insertion and removal

Available in three sizes: 6FR, 10FR and 14FR

Ordering Information

InterGuide



8070010

Code	Size	Diameter (OD)	Length	Box qty.
8070015	15FR	5.0mm	700mm	10
8070010	10FR	3.3mm	700mm	10
8070006	6FR	2.0mm	530mm	10

InterForm



8080010

Code	Size	Diameter (OD)	Length	Box qty.
8080014	14FR	4.7mm	340mm	10
8080010	10FR	3.3mm	340mm	10
8080006	6FR	2.0mm	275mm	10

IS6.5 • Issue 1 09.12