




Cirrus™ 2

A range of nebulisers and kits

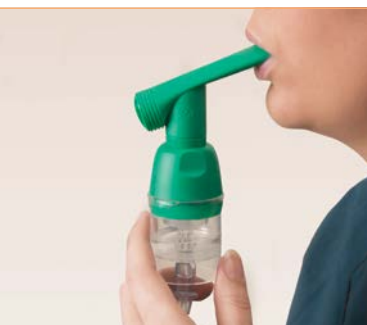

Oxygen and Aerosol Therapy • Nebuliser Therapy



Cirrus™2 nebuliser and nebuliser mask kits

The Cirrus2 nebuliser incorporates a number of features, including a graduated drug cup which provides guidance on how much drug is in the nebuliser. The cup has a quick-release fitting which enables the nebuliser to be attached, as well as detached, with a quarter turn. This speeds up the refilling process and provides a more reliable seal between the cup and the nebuliser body.

The Cirrus2 nebuliser is designed to deliver drugs for tracheobronchial deposition. At a driving gas flow of 8 L/min, the Mass Median Diameter (MMD) is 3.1 microns.

▶ Video available at www.intersurgical.com



Self-sealing T-piece

The Cirrus2 is available with a self-sealing T-piece, which enables the nebuliser to be attached (or detached) from the ventilator without interfering with patient ventilation or inadvertently activating any alarms.

▶ Video available at www.intersurgical.com





Cirrus™2 nebuliser, Intersurgical EcoLite™ mask kit

Welcome to the Comfort Zone

The unique design and materials used to manufacture the Intersurgical EcoLite mask provides improved patient comfort and a reduced environmental impact.

Comfortable for the Patient

Improved patient comfort has been key to the development of the Intersurgical EcoLite mask. The latest manufacturing technology has enabled us to combine two non-PVC materials in the same product, the polypropylene material forming the body of the Intersurgical EcoLite mask is clear, lightweight and rigid enough to maintain the masks shape, whilst the second softer TPE material is utilised in the seal, which is in contact with the patients face.

Comfortable for the Environment

The Intersurgical EcoLite mask is an important part of our Eco range, designed as part of our ongoing focus on sustainable development. Various groups have questioned the use of PVC in medical products and its impact on the environment. The utilisation of alternative materials has resulted in the elimination of PVC from the Intersurgical EcoLite mask, reducing the environmental impact.



More about our Eco range

www.intersurgical.com/info/eco

Video available at www.intersurgical.com



Incurved nose seal

Conforms to different nose shapes designed to prevent aerosolised drug entering patient's eyes

No metal nose clip

MRI compatible

Soft face seals

Increased patient comfort


Elastic can be positioned under or over the ears

Below ear position eliminates trauma to top of ears



A choice of "under chin" or "on chin" positions

Provides a better fit on a wider range of patient face shapes

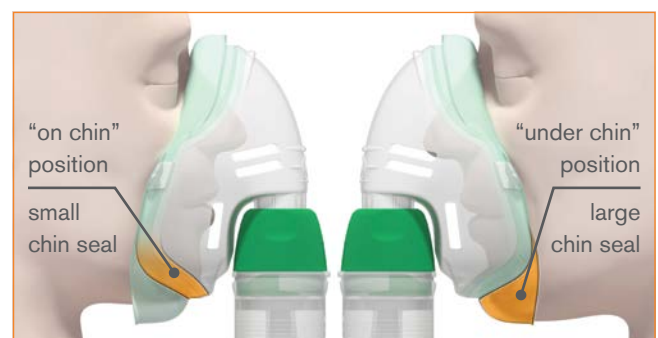



illustrated

Cirrus™2 nebuliser and nebuliser mask kits

Code	Description	Tube length	Box Qty.
1402000	Cirrus2 nebuliser		75
1453015	Cirrus2 nebuliser, adult, Intersurgical EcoLite mask kit with tube	2.1m	30
1453017	Cirrus2 nebuliser, adult, Intersurgical EcoLite mask kit with non-PVC tube	2.1m	30
1454015	Cirrus2 nebuliser, paediatric, Intersurgical EcoLite mask kit with tube	2.1m	36
1455000	Cirrus2 nebuliser, mouthpiece kit with tube, single patient use	2.1m	35
1464000	Cirrus2 nebuliser, universal mouthpiece T-kit with tube	1.8m	40
1464004	Cirrus2 nebuliser universal mouthpiece T-kit, filter and Sure Loc tube	1.8m	20
1464021	Cirrus2 nebuliser, mouthpiece kit, filter and tube	2.1m	20
2605000	Cirrus2 nebuliser breathing system T-kit, 22mm and tube	1.8m	40
2606000	Cirrus2 nebuliser breathing system T-kit, 15mm and tube	1.8m	35
2608000	Cirrus2 nebuliser breathing system T-kit, 10mm and tube	1.8m	45



illustrated

Self-sealing T-piece

Code	Description	Box Qty.
2609000	Cirrus2 nebuliser self-sealing T-piece kit, 22mm	35
1814000	Self-sealing in-line T-piece 22M-22F, 22M	40

Make an enquiry

Lower environmental impact product



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK
 T: +44 (0)118 965 6300 F: +44 (0)118 965 6356 info@intersurgical.com www.intersurgical.com



The manufacturer Intersurgical Ltd is certified to ISO 9001:2015, ISO 13485:2016 and ISO 14001:2015

Please think before you print
 Save energy and paper.
 If you must print this information sheet please print it double sided.

IS10.23 • Issue 12 07.19