



Intersurgical EcoLite™ adult high concentration oxygen mask

The lightest and most comfortable mask available



Oxygen and Aerosol Therapy ▪ Oxygen Therapy



Intersurgical EcoLite™ adult high concentration oxygen mask

Welcome to the Comfort Zone

The unique design and materials used to manufacture the Intersurgical EcoLite mask provides improved patient comfort and a reduced environmental impact.

Comfortable for the Patient

Improved patient comfort has been key to the development of the Intersurgical EcoLite mask. The latest manufacturing technology has enabled us to combine two non-PVC materials in the same product, the polypropylene material forming the body of the Intersurgical EcoLite mask is clear, lightweight and rigid enough to maintain the masks shape, whilst the second softer TPE material is utilised in the seal, which is in contact with the patients face.

▶ Video available at www.intersurgical.com

Incurved nose seal

Conforms to different nose shapes designed to prevent oxygen entering patient's eyes

No metal nose clip

MRI compatible

Two integrated safety valves

Allows patient to breathe room air in the event oxygen supply failure

A choice of "under chin" or "on chin" positions

Provides a better fit on a wider range of patient face shapes




Comfortable for the Environment

The Intersurgical EcoLite mask is an important part of our Eco range, designed as part of our ongoing focus on sustainable development. Various groups have questioned the use of PVC in medical products and its impact on the environment. The utilisation of alternative materials has resulted in the elimination of PVC from the Intersurgical EcoLite mask, reducing the environmental impact.



More about our Eco range

www.intersurgical.com/info/eco

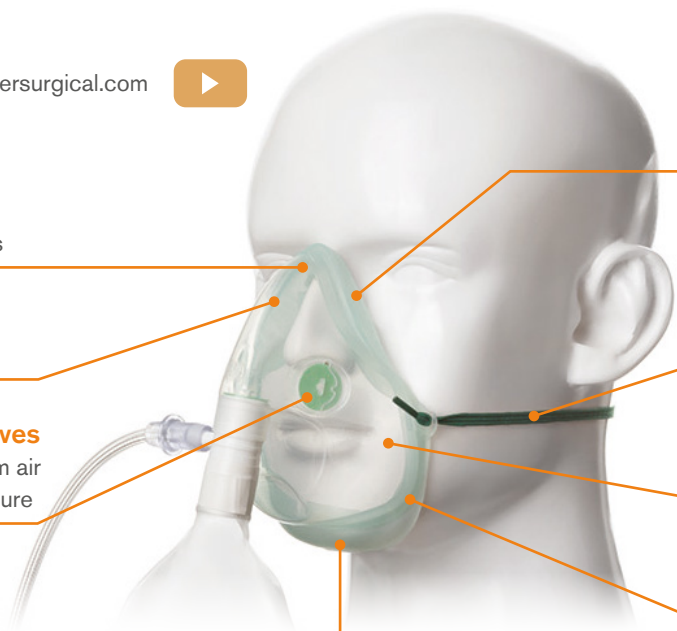


Soft face seals
Increased patient comfort

Elastic can be positioned under or over the ears
Below ear position eliminates trauma to top of ears

Clear mask
permits visualisation of the patient's nose and lips

Improved sealing
Higher oxygen concentrations

Oxygen percentage delivered Intersurgical EcoLite™ compared to PVC high concentration masks

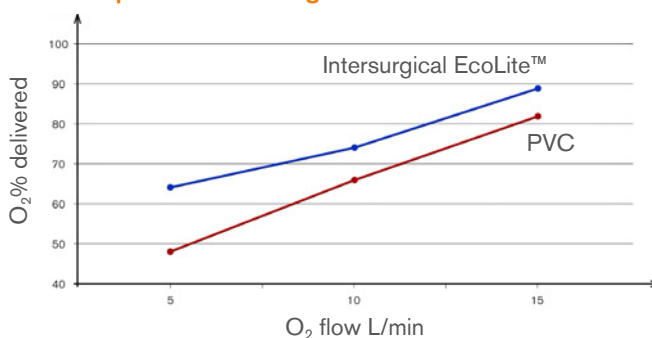

Code	Description	Tube length	Box Qty.
1181015	Intersurgical EcoLite, adult, high concentration oxygen mask with tube	2.1m	24

Lower environmental impact product [Make an enquiry](#)

INTERSURGICAL COMPLETE RESPIRATORY SYSTEMS Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK
T: +44 (0)118 965 6300 F: +44 (0)118 965 6356 info@intersurgical.com www.intersurgical.com

The manufacturer Intersurgical Ltd is certified to ISO 9001:2015, ISO 13485:2016 and ISO 14001:2015 **Please think before you print** Save energy and paper. If you must print this information sheet please print it double sided. IS10.17 • Issue 7 07.19